
SANTÉ

Bruksanvisning

COMFORT Go
Massagekudde

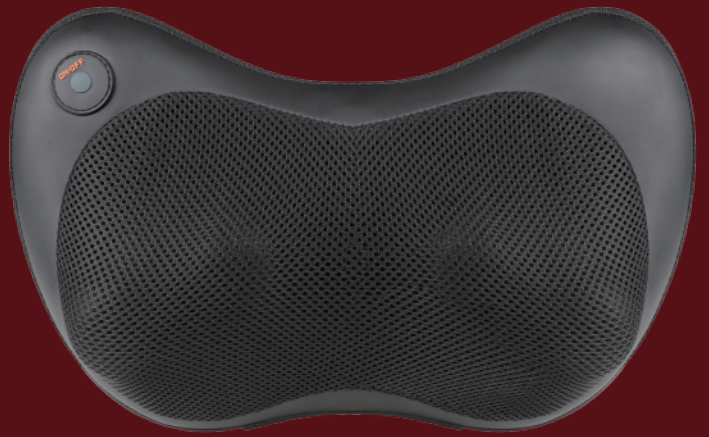SANTÉ

Santé Sweden AB
Telefon: +46 20 27 27 20
www.massageprodukter.se

Vi gratulerar till ditt köp av Massagekudde COMFORT Go.

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
Instruction manual

COMFORT Go
Massage pillow

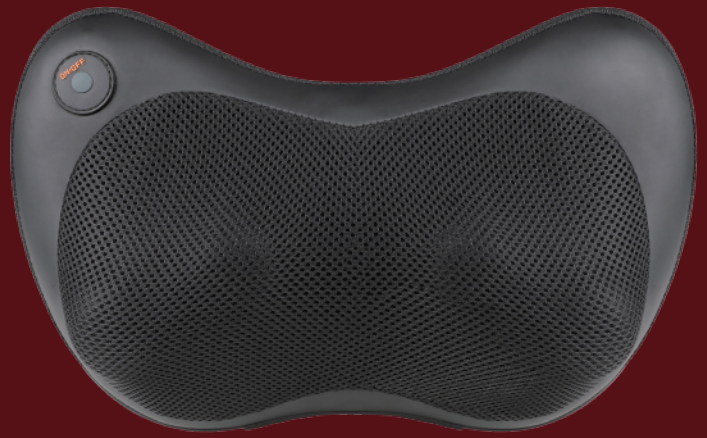Congratulations on your purchase of Maassage pillow COMFORT Go.

By choosing this massage product from Santé, you have taken a conscious step toward better health, enhanced well-being, and a more balanced everyday life. We are proud to be part of your journey toward recovery and improved quality of life.

Santé is committed to providing innovative and high-quality products in the field of health and wellness. With a focus on modern technology, thoughtful design, and user-friendliness, we develop solutions that help you care for both body and mind—whether for relaxation after a long day or as part of your regular self-care routine.

Our products are developed with care and carefully selected components to ensure functionality, durability, and comfort. We continuously strive to improve and further develop our range to meet our customers' needs and expectations.

To ensure you get the most out of your product and can enjoy its benefits for a long time, we recommend that you read and follow the user manual carefully. This will help ensure optimal use and a safe experience.

Thank you for choosing Santé — we wish you many relaxing and rejuvenating moments!

CONTENTS

Description ----------/page 1

Safety precautions ----------/page 1

Product illustration----------/page 2

Operation instructions----------/page 3

Application ----------/page 3

Benefits ----------/page 3

Technical data ----------/page 4

Description

Whether you're looking for a quick relaxing fix from the hustle and bustle of everyday life or some serious relief from the stubborn knots that just won't quit on your back, neck, etc. the COMFORT Go is the electric massager that banishes the aches for relaxation anytime, anywhere. This portable massager with heat therapy makes for excellent relief on your neck, shoulders, back, legs, arms, or thighs, for that full-body massage therapy in the comfort of your own home.

Safety Precautions

Important safety Warning

- Do not insert debris into the gap between the massage head and the shell.
- Do not splash water on the product to avoid accidents caused by short circuit.
- If the power cord or plug is damaged, do not disassemble by yourself and do not repair it without authorization.
- Please contact the dealer in time for maintenance.
- Please use a power source matching this appliance.
- Please unplug the power plug after use or before cleaning to avoid causing equipment or personal damage.
- Please follow the instructions in this manual to use this massager.
- Do not use unrecommended accessories and accessories.
- Do not use this massager outdoors.
- Do not use the appliance in case of damaged enclosure or power cord.
- Do not use this massager when you are drunk or feel unwell.
- Do not use the appliance with excessive force to avoid injury.
- Do not use this product within one hour after meals.
- For patients with hand injury or recent hand operation, please follow the doctor's instructions and use carefully.


Use environment

- Do not use the massager in a highly damp environment such as the bathroom.
- Do not use the massager immediately when the ambient temperature changes sharply.
- Do not use the massager in severe flying dust or corrosive gas environment.

People who are not suitable to use this device

- People with skin infections.
- Osteoporosis patients.
- Patients with heart problems and wearing electronic medical devices such as pacemakers.
- People with high fever.
- Pregnant or menstruating women.
- People with physical trauma or surface diseases.
- Children under the age of 14 and those who are insane are not allowed to use this product without supervision.
- Require rest or feel unwell as advised by the doctor.
- Do not use this massager when the body is wet.
- The appliance has a heating surface and please attention and carefully use to those people who are not sensitive to heat.

Product illustration

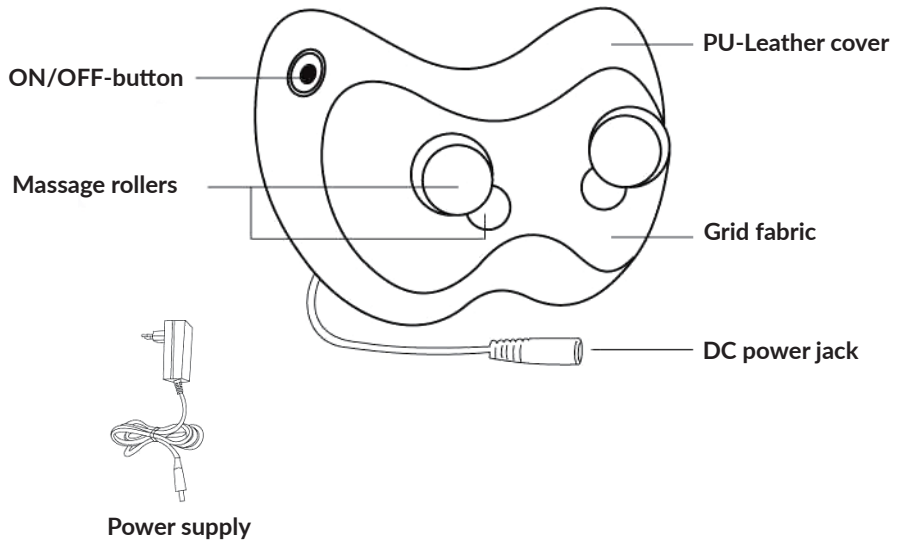Operation instructions

1. Make sure the massager is operated under safe environment.
2. Insert the power plug into the car plug adapter or wall socket.
3. Press the ON/OFF-button to start the massager and thermal therapy function.
4. The massage rollers will rotate in clockwise and counterclockwise, and the rotation direction will change automatically every 1 minute.
5. Double press the ON/OFF-button, and the massage rollers will change the rotation direction and the manual steering is on.
6. Press the ON/OFF-button to stop the thermal therapy function and press the button again to start the function.
7. Press and hold the ON/OFF-button for 3s to turn off the massager.
8. Unplug the massager.

Application

- Neck
- Back
- Stomach
- Thighs
- Arms
- Calves

Benefits

- Improve blood circulation
- Alleviate muscle tension/aches
- Improve metabolic system
- Relieve stress and fatigue

Technical data:
Product type: Massage pillow
Product name: COMFORT Go
Article no. MP-P11C-N-1
Output voltage: DC 9 V
Power consumption: 18 W
Timer: 15 min.

Santé Sweden AB
VAT: SE556646715401
Phone: +46 20 27 27 20
www.massageprodukter.se

EU Declaration of Conformity

Compliant with European regulatory standards EMC 2014/30/EU and LVD 2014/35/EU.

Warranty Conditions:

The retailer guarantees a statutory period of warranty for the product described on the reverse. Warranty is valid from the date of purchase. The purchase date must be documented with proof of purchase. The retailer will repair or replace free-of-charge all items covered in each market for which faults have been discovered. A prerequisite is that you send the faulty device with proof of purchase before the guarantee or warranty period expires. The warranty claim shall not be applicable if a fault is established which has arisen e.g. due to external influences or as a consequence of a repair or modification which was not carried out by the manufacturer or an authorised dealer.

The following points are not part of the warranty:

- Wear parts, e.g. abrasives, working synthetical parts, fabric covers, overlay / cover
- Wear caused by improper use
- Commercial use e.g. exhibition use, repeated use without intermediate breaks of 40 minutes
- Exposing the device to heat, humidity of liquids
- Damage due to sharp and/or sharp-edged objects
- Wanton destruction, misappropriation, overloading or accidental damage
- Damage and soiling by domestic animals
- Contamination caused by body or hair care products
- Damage caused by use of acids, solvent or alcohol agents

The warranty granted by the retailer is limited to repair and exchange of the product. Under this warranty, the manufacturer or the retailer shall have no further liability and are not responsible for damages which have arisen due to disregard of the operating instructions and /or improper use of the product. Where upon examination of the product by the manufacturer it is found that the warranty request does not concern a defect covered by the warranty or that the warranty period has expired, the customer must bear the costs for the inspection and repair.As manufacturer, Casada ensures the delivery of faultless goods to buyers/partners. Proven production faults are entitled to an exchange after registration and confirmation.

Legal note: All contents do not contain healing statements. The diagnosis and treatment of illnesses and other physical disorders require a treatment by a physician, non-medical practitioner or therapist. The statements are strictly informative and may not be used as a substitute for a professional medi-cal treatment. Every user is beholden to ascertain by accurate examination and, if required, after consultation of a specialist, if the application is the right one for their specific cases. Each dose or application takes place at the risk of the user. We dissociate ourselves from any healing statements or promises.

Copyright pictures and text are subject to the copyright of Santé Sweden AB. and may not be used without explicit confirmation. Copyright © 2026 Santé Sweden AB. All rights reserved.

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Upphovsrättsskyddade bilder och text omfattas av Santé Sweden AB:s upphovsrätt och får inte användas utan uttrycklig bekräftelse. Copyright © 2026 Santé Sweden AB. Alla rättigheter förbehålls.