SANTÉ

Bruksanvisning

PULSE Mini
Massagepistol

SANTÉ

Santé Sweden AB
Telefon: +46 20 27 27 20
www.massageprodukter.se

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Disposal and Recovery
This product should not be discarded with household waste. When you decide to dispose of this product, please follow the local environmental laws and guidelines.

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SANTÉ

Instruction manual

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Congratulations on your purchase of Massage gun PULSE Mini.

By choosing this massage product from Santé, you have taken a conscious step toward better health, enhanced well-being, and a more balanced everyday life. We are proud to be part of your journey toward recovery and improved quality of life.
Santé is committed to providing innovative and high-quality products in the field of health and wellness. With a focus on modern technology, thoughtful design, and user-friendliness, we develop solutions that help you care for both body and mind—whether for relaxation after a long day or as part of your regular self-care routine.
Our products are developed with care and carefully selected components to ensure functionality, durability, and comfort. We continuously strive to improve and further develop our range to meet our customers' needs and expectations.
To ensure you get the most out of your product and can enjoy its benefits for a long time, we recommend that you read and follow the user manual carefully. This will help ensure optimal use and a safe experience.
Thank you for choosing Santé — we wish you many relaxing and rejuvenating moments!

CONTENTS

Attention----------/page 2
Please read before operation----------/page 3
Product spercification----------/page 4
Product feature----------/page 4
How to operate ----------/page 5
Faulty inspect and reset----------/page 8
Battery dismantling steps----------/page 9
Storage and maintenance----------/page 10
Product information----------/page 11

ATTENTION

In order to operate well and safely , please read the Instructions carefully before use.
The user shall be responsible for any consequences caused by improper use or non-compliance with the instructions

1. Do not use the devise if you are drunk.
2. Please remove ornaments or objects that may affect when you massaging.
3. Do not use in the bathroom or other humid environment.
4. Keep the product away from water , fire and corrosive environment.
5. Please comply with the safety regulations when using in the specific environment.
6. Please do not use on the surgical site of your body.
7. Do not allow to use on those people that feel skin damage, redness, inflammation, congestion and other uncomfortable symptoms.
8. Do not allow to use on those patients with heart disease or recorded, acute diseases, skin diseases, malignant tumors, or any other unhealthiness.
9. Do not allow to use on those patients with blood diseases such as bleeding tendency, osteoporosis, soft tissue diseases or stroke.
10. Please do not use if you have slow movement, physical and mental dysfunction or intellectual disability. please follow the doctor's advice when you need this devise.
11. Please stop using if you feel uncomfortable during use.
12. Do not use massagers during pregnancy or period.
13. It is prohibited for children under 6 years old , children 6-12 years old must be used under the guidance and care of parents or persons with supervisory ability.
14. If you need to clean the appearance, use alcohol cotton to clean,No water.
15. Do not use sharp objects to scratch or stab the product.
16. As household electronic healthcare product, not as a medical device, and not suitable for diagnosis, treatment and other commercial purposes.
17. Please use the charging equipment offer by the seller to charge the product, and please be careful from electricity.
18. If there is any abnormal phenomenon during use, please cut off the power immediately. Do not disassemble or repair by yourself.
19. Some of the machines have the function of hot moxibustion, and the instruments would be heat on surface. People who are not sensitive to heat must be pay attention when using them

ATTENTION

- this appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved

- children shall not play with the appliance

- cleaning and user maintenance shall not be made by children without supervision

PLEASE READ BEFORE OPERATION

You're supposed to know
※ About frequency of use
Use it once a day in the morning and in the evening, and keep the time within the rated time.


※ Pay attention before use
1. Do not operate the product with your wet hands.
2. Please remove the wearing jewelry to avoid discomfort during massage.
3. Discontinue use and consult a physician if you feel pain or heat.

※ What will happen after use
The massaged area may turn red and leave a mark for a short time, but it will recover after a period of time.

※ Others would happen
During use, there will be a "buzz", "pop" and "boom" sound. This is the sound of adjustment, not an abnormality.

PRODUCT FEATURES

Accessories of this product include PULSE mini massage gun, four massage heads, charging cable, instruction.

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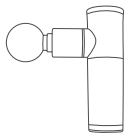Four massage heads

Charging cable

User Manual

PRODUCT FEATURE

Compact and portable

HOW TO OPERATE

Preparation before use

1. Check whether the products and accessories are complete. Do not use them if they are damaged.

2. Check whether the battery is sufficient. If the battery is low, please charge it before using it.

Charging method

1. Insert the Type-C plug on the USB cable into the Type-C socket of the base, and connect the other end to the adapter USB interface.

2. Put the host machine on the base.

The host indicator light is always on pink, indicating charging is in progress, and the green light is always on, indicating charging is completed.

※ It is normal for the devise and power adapter to heat up slightly during charging.

Charging time

1. Charging time is 4-5 hours.

2. Use 8 times for 10 minutes at a time when fully charged.

(Due to the different massage intensity, there will be differences)


3. When the indicator light flashes red, please charge it immediately.

Electricity

The machine will shut down automatically in low power.


When the charging is over

Unplug the power adapter from the devise and socket when finish charging.

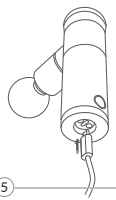HOW TO OPERATE


Schematic diagram of product function configuration

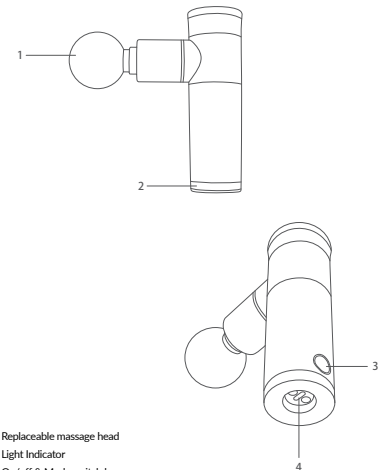1. Replaceable massage head
2. Light Indicator
3. On/off & Mode switch key
4. Charging interface

HOW TO OPERATE

Operational approach

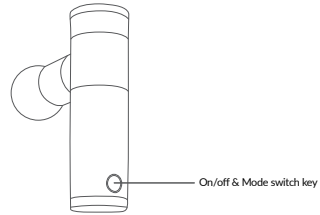This product cannot be used while charging.Please unplug the power adapter and use it after charging.

1. Install or replace the massage head.Press the power key.

2. On/off machine: Long press the "on/off" button for two seconds, the machine will start up, the indicator light will be light green, and the machine will run at a slow speed by default.

3. Intensitive adjustment: continuousing short press the on/off key to adjust five speed from weak to strength and change from strong to weak cycle adjustment , different color to indicate different intensitive: the first level is green light, a second level emerald, three level is purple, four level "light pink, firth level sky blue.

4. Aim to the parts that need to be massaged.
※ Excessive force will cause the self-protection and automatically stop working. After reaching 3 times of self-protection, the device will turn off. Then restarting, the device will work normally.

5. Automatic shutdown: After 10 minutes of use, the product will automatically shut down.

6. Shutdow by hand : Long press the "on/off" button for two seconds.

FAULTY INSPECT AND RESET

Please check if there is a fault in the operation. If the fault cannot be removed or the product still cannot be operated satisfactorily, please cut off the power supply and be repaired by professional technicians. Please do not disassemble the product by yourself!

Basic faults	Possible causes	Basic troubleshooting methods
Unable to boot	**is the battery low**	**if yes, please charge it**
Weak massage	**is the battery low**	**if it is, please charge it immediately**

BATTERY DISMANTLING STEPS

The lithium ion battery used in this product is a valuable resource that can be recycled and regenerated. Please take out the lithium ion battery after use and give it to a professional organization that can recycle and regenerate rechargeable batteries for processing.

Steps to remove the battery:

1. Cut off the external power supply and remove the product shell;
2. Unplug batteries;
3. Remove the battery;
4. Waste batteries should be collected tothe recycling station;
5.Please follow the local environmental protection laws and discard guidelines and protect the environment.

Note:

- The battery must be removed from the appliance before it is discarded;
- The appliance must be powered off when the battery is removed;
- Batteries should be disposed of safely.

⚠ Non-professional personnel please do not disassemble by yourself!

STORAGE AND MAINTENANCE

Storage

storage over than 3 months, the machine should be recharged after unpacking and inspection. Do not place it close to the floor, walls or roof. Indoor should be well ventilated, avoid intense sunlight illuminate and corrosive gas erosion.

Transport

The product package is equipped with simple anti-shake, suitable for air, railway, road and ship transportation, should avoid rain, snow splashing, inversion and collision.

Power inspection

Check the power supply. If the power supply voltage exceeds the specified range, do not start up and use.

Maintenance

Take measures to prevent moisture, rodents, and pollution. The cleaning should be done in the shutdown state. After use, place the instrument in a safe place, and never put heavy objects on the instrument. When a fault occurs, please ask professional and technical personnel to repair, non-professionals are not allowed to disassemble the product by themselves.

Disposal and Recovery

This product should not be discarded with household waste. When you decide to dispose of this product, please follow the local environmental laws and guidelines.

PRODUCT INFORMATION

Technical data:

Product type: Massage gun
Product name: PULSE Mini
Article no: MP-A28-N-1
Rated voltage: 5 V ⎓
Rated power: 10 W
Battery capacity: 2200mAh
Product weight: 0.5 kg
Measures: L 47 x W 90 x H 152mm

Santé Sweden AB
VAT: SE556646715401
Phone: +46 20 27 27 20
www.massageprodukter.se

EU Declaration of Conformity

Compliant with European regulatory standards EMC 2014/30/EU and LVD 2014/35/EU.

Warranty Conditions:

The retailer guarantees a statutory period of warranty for the product described on the reverse. Warranty is valid from the date of purchase. The purchase date must be documented with proof of purchase. The retailer will repair or replace free-of-charge all items covered in each market for which faults have been discovered. A prerequisite is that you send the faulty device with proof of purchase before the guarantee or warranty period expires. The warranty claim shall not be applicable if a fault is established which has arisen e.g. due to external influences or as a consequence of a repair or modification which was not carried out by the manufacturer or an authorised dealer. The following points are not part of the warranty:

- Wear parts, e.g. abrasives, working synthetical parts, fabric covers, overlay / cover
- Wear caused by improper use
- Commercial use e.g. exhibition use, repeated use without intermediate breaks of 40 minutes
- Exposing the device to heat, humidity of liquids
- Damage due to sharp and/or sharp-edged objects
- Wanton destruction, misappropriation, overloading or accidental damage
- Damage and soiling by domestic animals
- Contamination caused by body or hair care products
- Damage caused by use of acids, solvent or alcohol agentsThe warranty granted by the retailer is limited to repair and exchange of the product. Under this warranty, the manufacturer or the retailer shall have no further liability and are not responsible for damages which have arisen due to disregard of the operating instructions and /or improper use of the product. Where upon examination of the product by the manufacturer it is found that the warranty request does not concern a defect covered by the warranty or that the warranty period has expired, the customer must bear the costs for the inspection and repair.As manufacturer, Casada ensures the delivery of faultless goods to buyers/partners. Proven production faults are entitled to an exchange after registration and confirmation.

Legal note: All contents do not contain healing statements. The diagnosis and treatment of illnesses and other physical disorders require a treatment by a physician, non-medical practitioner or therapist. The statements are strictly informative and may not be used as a substitute for a professional medi-cal treatment. Every user is beholden to ascertain by accurate examination and, if required, after consultation of a specialist, if the application is the right one for their specific cases. Each dose or application takes place at the risk of the user. We dissociate ourselves from any healing statements or promises.

Copyright pictures and text are subject to the copyright of Santé Sweden AB. and may not be used without explicit confirmation. Copyright © 2026 Santé Sweden AB. All rights reserved.

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Upphovsrättsskyddade bilder och text omfattas av Santé Sweden AB:s upphovsrätt och får inte användas utan uttrycklig bekräftelse. Copyright © 2026 Santé Sweden AB. Alla rättigheter förbehålls.