SANTÉ

Bruksanvisning

PULSE Pro
Massagepistol

SANTÉ

Santé Sweden AB
Telefon: +46 20 27 27 20
www.massageprodukter.se

Vi gratulerar till ditt köp av Massagepistol PULSE Pro.

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Instruction manual

PULSE Pro
Massage gun

Congratulations on your purchase of Massage gun PULSE Pro.

By choosing this massage product from Santé, you have taken a conscious step toward better health, enhanced well-being, and a more balanced everyday life. We are proud to be part of your journey toward recovery and improved quality of life.

Santé is committed to providing innovative and high-quality products in the field of health and wellness. With a focus on modern technology, thoughtful design, and user-friendliness, we develop solutions that help you care for both body and mind—whether for relaxation after a long day or as part of your regular self-care routine.

Our products are developed with care and carefully selected components to ensure functionality, durability, and comfort. We continuously strive to improve and further develop our range to meet our customers' needs and expectations.

To ensure you get the most out of your product and can enjoy its benefits for a long time, we recommend that you read and follow the user manual carefully. This will help ensure optimal use and a safe experience.

Thank you for choosing Santé — we wish you many relaxing and rejuvenating moments!

CONTENTS

Important safeguards ----------/page 2-3

Product description ----------/page 4-5

Charging ----------/page 6

Accessories listing ----------/page 7

Operating instruction ----------/page 8

Maintenance, cleaning & storage ----------/page 8-09

Product parameters ----------/page 9

Usage ----------/page 10-15

After-sales service ----------/page 16-18

Fault exclude ----------/page 19

Important Safeguards

⚠ Statement
To reduce the risk of malfunction, electric shock, injury, fire, death, and equipment damage, please read all instructions and warnings on this product in this manual before using this product.

Prohibited

1. It is forbidden for minors to use. If injured, it is forbidden to use this equipment and please consult your doctor before using this device.It is forbidden to use on non-dry, clean skin surface. It can be used across the clothes and gently pressed and moved on the skin for about 60 seconds.
2. Do not use above the neck or any hard or boned part of the body, or sensitive areas. Use the product only on the soft tissues of the body without causing pain or discomfort.
3. Do not massage areas that are bruised or on open wounds, stop use if you have any signs of abnormal pain or discomfort.
4. Do not place any objects on the vents of the product.
5. Do not immerse in water or allow water to enter the vents of the product.
6. Do not drop or close the product to fire source.
7. Do not drop or improperly use the product.
8. Do not use a non-DC 26 volt charger to charge this product.
9. Do not disassemble the screws or attempt to disassemble the
10. equipment and modify the fascia gun in any way.
11. Do not allow your fingers, hair, or other parts of your body to get close to the shaft and back of the vibrating head as they may get caught.

Do not touch this area when the device is running without massage head installed.

Do not touch this area when the device is running.

Warnings and Reminders

1. Keep away from liquids or heat sources. Keep the vents free of dust and debris.
2. After charging is complete, remove the charger from the battery.
3. Please mind that the device can not be used when charging.
4. Please use the device according to the instructions for use.
5. Carefully inspect the equipment before each use.
6. In order to ensure that the product can be used continuously, please charge the product in time to ensure that the product is always powered, and prevent the product from being unable to turn on due to lack of power.


Important

If you want to use it in any of the following situations, you should not use a fascial massage gun without your doctor's consent.


Pregnancy, complications of diabetes (such as neuropathy or retinal damage), wearing a pacemaker, recent surgery, epilepsy or migraine, disc herniation, advancement of the spine, vertebral detachment, or rigid joints of the spine, recent joints Replacement surgery or IUD, metal needle or steel plate any problems related to your health. These contraindications do not mean that you are not able to use vibration equipment or training equipment, we recommend that you consult a doctor before use.

The ongoing research is about the impact of vibration training on specific medical conditions. This is likely to reduce the types of contraindications mentionedabove. Practical experience has shown that in many cases, the integration of vibration training into a treatment regimen is also a therapeutic recommendation. This must be done with or with the advice of a doctor, expert or physical therapist.


Product Description

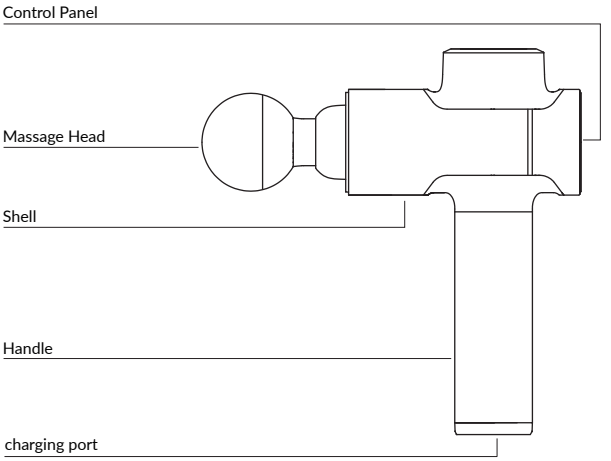1. High-speed and brush-less balance motor with dynamic balance mute technology.
2. Various massage heads aim at different areas of body, which also reduce damage to bones.
3. Rechargeable lithium-ion battery (4+ Hours per Charge).
4. Handle vertical design, grip more comfortable(48mm).
5. Five-speed settings adjust free.
6. User-friendly design of sweatproof.
7. Battery and body integrated design, keep handle without risks of looseness.

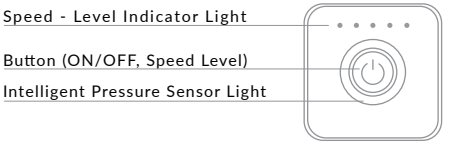1. Speed- Level Indicator Light:From Left to Right (1-5 level)
2. Battery -Level Indicator:4 Lights(100% full of power);
 3 Lights(65%-80% of power); 2 Lights(40%-60% of power);
 1 Light (Less power, please charge asap).
3. Intelligent pressure sensor light: As the resistance increases/decreases, The stronger the resistance, the fuller the lamp.

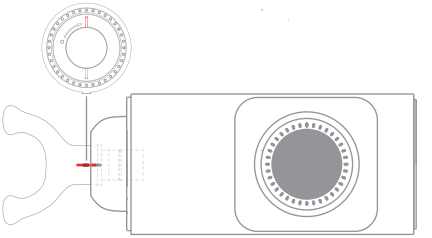How to install:

Install the massage head, rotate it clockwise (about 90°), until the vertical line is lined as the figure shown, indicating the massage head is in place.


Applicable age:

For adult use only. Pregnant woman should avoid use. See safety instructions.

Charging

1. **Fully Charge the battery for six hours before first use.**
2. **To charge, connect the DC end of the supplied 26 volt adapter to the battery charging port and plug the AC side cable into the wall outlet.**
3. **The Battery - Level Indicator is flashing regularly when the device is charging.**
4. **When the first light is on, the remaining power is 20%-35%; when the second light is on, the remaining power is 40%-60%; when the third light is on, the remaining power is 65%-80%; when the fourth light is on, the remaining power is 85%-100%.**
5. **The battery can be charged in the situation of any time and any quantity of electricity, the device can not be started when charging.**
6. **The battery is exhausted to less than 5% (the first left Battery - Level Indicator LED flashes) and needs to be recharged.**
7. **The average running time is about 2-3 hours depending on the speed level and pressure during use.**


Accessories listing

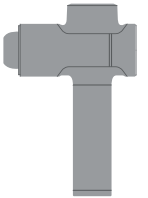Device

Carrying Case

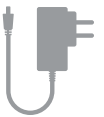Power Adaptor

User's Manual

Fork

Round

Flat

Bullet

Strawberry

Operating Instruction

1. Press the ON/OFF Switch to turn to the device to level 1. Each press of the button will increase one level until the level 5 speed. Touch the button at level 5 to enter the standby model (level 0), level 5 to level 0. Touch and hold again to turn the device to level 1, and so on.
2. Under any speed level, press and hold can close the device Breathing light flashing.
3. In the standby state, press and hold to directly shut down, 15 seconds without operation will automatically shut down.
4. The device will stop running after 10 minutes of continuous use for intelligent protection. Restart again if you use it again after the device shut down.
5. The massage head stop running after external force, This is a abnormal operation protection function setting, from running to standby mode, short press the button to restart, it's normal situation.

Maintenance, Cleaning, Storage & Transportation

1. Use a slightly damp towel to clean device and a dry with a soft cloth.
2. For storage or travel, store the device in supplied portable bag.
3. Store in a dry, well-ventilated with no corrosive gas place, no less than 250mm from the ground, and no more than 5 stacked layers.
4. If not been used for a long time, the device should be cleaned then packaged in storage, and charged every one month.
5. This product is suitable for automobile, train, ship and airplane transportation.
6. This product should be fixed in transportation to prevent violent vibration, collision, tumbling and falling, and should be put upward, put down gently and avoid wet or damp areas.

Environmental Protection

At the end of the service life of this product, please do not throw it away with general household waste. Put this product in the designated recycling station of the government. You can also send this product back to the after-sales service center of the factory. The personnel of the service center will handle it with environmentally friendly way.

Technical data:
Product type: Massage gun
Product name: PULSE Pro
Article no. MP-A59B-N-1
Power supply: 100–240 V ~ 50/60 Hz.
Interface: DC 26 V / 0.9 A
Battery: Rechargeable Litium-Ion battery, 2600 mAh.
Charging time: Approx. 4 hours
Speed levels: 5 adjustable levels
Weight: 1,2 kg.
Measures: 164 × 71 × 238 mm.
Useage environment: 0 – 40 °C

Santé Sweden AB
VAT: SE556646715401
Phone: +46 20 27 27 20
www.massageprodukter.se

Usage

Fork

For lumbar vertebrae and large muscle groups

Round

For large muscles groups

Flat

For various parts of the body


Bullet

For deep muscle groups


Strawberry

For sensitive parts of muscle or older person

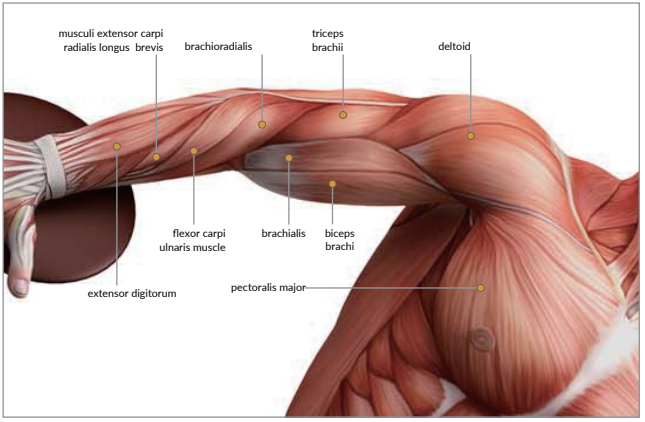Arm


brachioradialis	10-15s x3 times
flexor carpi radialis	10-15s x3 times
palmaris longusmuscle	12-15s x3 times
flexor carpi ulnaris muscle	10-15s x3 times
extensor digitorum	05-10s x2 times
musculi extensor carpi radialis longus	05-10s x2 times
musculi extensor carpi radialis brevis	05-10s x2 times
biceps brachii	12-20s x3 times
brachialis	12-20s x3 times
deltoid	12-20s x3 times
pectoralis major	12-20s x3 times

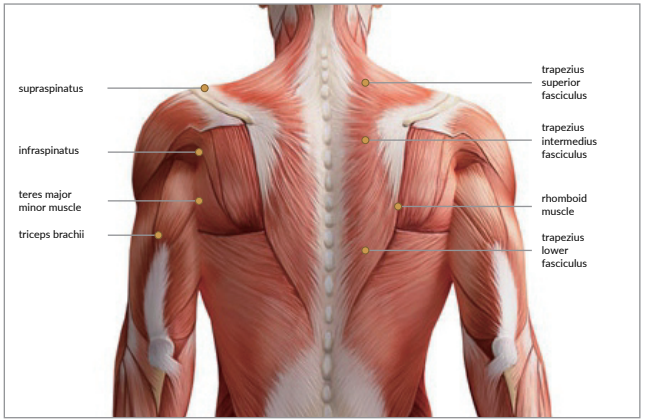Shoulder

triceps brachii	15-20s	x3 times
rhomboid muscle	20-30s	x3 times
trapezius superior fasciculus	15-20s	x3 times
trapezius intermedius fasciculus	15-20s	x3 times
trapezius lower fasciculus	15-20s	x2 times
supraspinatus	10-15s	x2 times
levator scapulae muscle	10-12s	x2 times
teres major muscle	10-12s	x3 times
teres minor muscle	10-12s	x3 times
infraspinatus	10-12s	x3 times


Tips:Massage the muscles from middle to sides, from top to bottom.

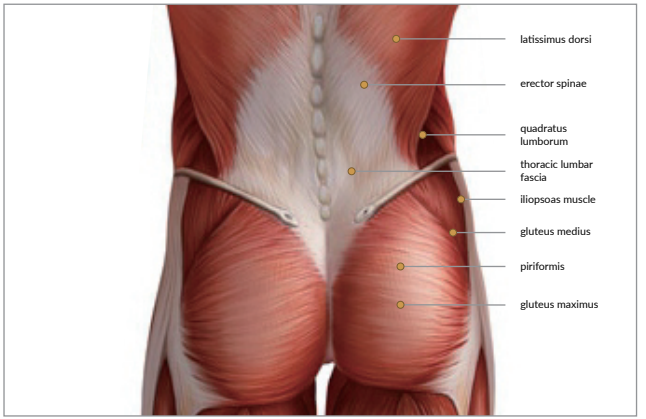Waist

erector spinae（lumbar segments）	30-60s	x3 times
latissimus dorsi	30-60s	x3 times
thoracic lumbar fascia	40-60s	x3 times
psoas major	30-60s	x3 times
internal and external oblique muscle	10-15s	x2 times
erector spinae（thoracic segments）	15-20s	x2 times
quadratus lumborum	20-30s	x2 times
iliopsoas muscle	20-30s	x2 times
gluteus maximus	30-60s	x3 times
gluteus medius	30-60s	x3 times
piriformis	60-90s	x3 times


Tips: Massage the muscles from middle to sides, from top to bottom.
Can follow the muscle direction on large muscle groups.

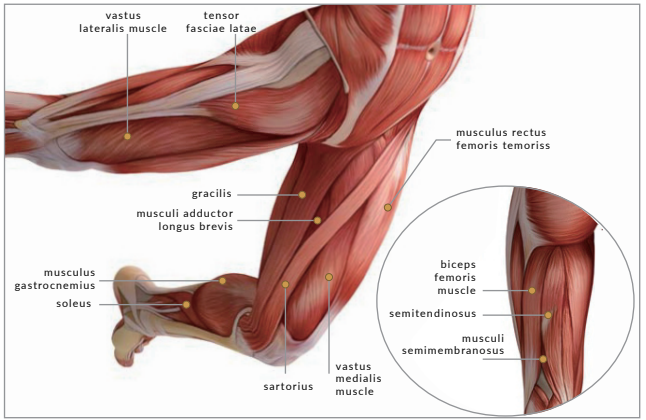Leg

biceps femoris muscle	20-30s x3 times
semitendinosusi	20-30s x3 times
musculi semimembranosus	20-30s x3 times
gracilis	10-15s x3 times
musculi adductor brevie	10-15s x2 times
musculi adductor longus	10-15s x2 times
musculus rectus femoris temoriss	30-45s x2 times
vastus lateralis muscle	30-45s x2 times
vastus medialis muscle	30-45s x3 times
sartorius	30-45s x3 times
sartorius	15-20s x3 times

Tips: Massage the muscles from top to bottom.

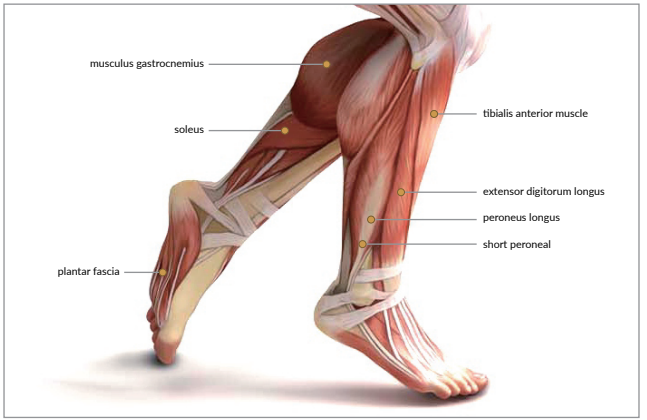Foot

musculus gastrocnemius	20-30s x3 times
soleus	20-30s x3 times
tibialis anterior muscle	15-20s x3 times
peroneus longus	15-20s x3 times
short peroneal	15-20s x2 times
plantar fascia	20-30s x2 times

Tips: Massage the muscles from top to bottom.

User info

Name: _____ **Purchasing date:** _____

E-mail: _____ **Contact number:** _____

Address: _____

Product information

Product model: _____ **Identification of product:** _____

Fault cause: _____

Dealer information

Dealer name: _____ **Contact number:** _____

Invoice No: _____ **Sales date:** _____

Address: _____

Massage Gun Operational Troubleshooting

Fault Description	Fault Exclude
1.Massage head slides out and drops	Please check that the massage head is installed in place (See instruction for detailed operation in 5th page)
2.Device unable to start	Please check that if the battery is adequate or not(Please charge over 8 hours if the device is not been used for a long time)
3.Device unable to charge	1. Please check whether the charging port and power supply are plugged in 2. Check whether the voltage and current of the charging head are matched with the product requirements
4.Device suddenly stop after running for 10 minutes	This is a long time running protection function setting, from running to standby mode, short press the button to restart, it's normal situation. (See instruction for detailed operation in 8th page)
5.The massage head stop running after external force	This is a abnormal operation protection function s etting, from running to standby mode, short press t he button to restart, it's normal situation (See instruction for detailed operation in 8th page)

Reasonable use of the product can effectively extend the service life, Suggestion:

1. Do not throw or drop when using.
2. The maximum continuous using time is not more than 30 minutes,The maximum continuous use time should not exceed 30 minutes. Keep away from damp environment when using or storing.
3. If you don't use it for a long time, please charge the product once every 1 month to ensure that the product is always powered.

EU Declaration of Conformity

Compliant with European regulatory standards EMC 2014/30/EU and LVD 2014/35/EU.

Warranty Conditions:

The retailer guarantees a statutory period of warranty for the product described on the reverse. Warranty is valid from the date of purchase. The purchase date must be documented with proof of purchase. The retailer will repair or replace free-of-charge all items covered in each market for which faults have been discovered. A prerequisite is that you send the faulty device with proof of purchase before the guarantee or warranty period expires. The warranty claim shall not be applicable if a fault is established which has arisen e.g. due to external influences or as a consequence of a repair or modification which was not carried out by the manufacturer or an authorised dealer. The following points are not part of the warranty:

- Wear parts, e.g. abrasives, working synthetical parts, fabric covers, overlay / cover
- Wear caused by improper use
- Commercial use e.g. exhibition use, repeated use without intermediate breaks of 40 minutes
- Exposing the device to heat, humidity of liquids
- Damage due to sharp and/or sharp-edged objects
- Wanton destruction, misappropriation, overloading or accidental damage
- Damage and soiling by domestic animals
- Contamination caused by body or hair care products
- Damage caused by use of acids, solvent or alcohol agents

The warranty granted by the retailer is limited to repair and exchange of the product. Under this warranty, the manufacturer or the retailer shall have no further liability and are not responsible for damages which have arisen due to disregard of the operating instructions and /or improper use of the product. Where upon examination of the product by the manufacturer it is found that the warranty request does not concern a defect covered by the warranty or that the warranty period has expired, the customer must bear the costs for the inspection and repair.As manufacturer, Casada ensures the delivery of faultless goods to buyers/partners. Proven production faults are entitled to an exchange after registration and confirmation.

Legal note: All contents do not contain healing statements. The diagnosis and treatment of illnesses and other physical disorders require a treatment by a physician, non-medical practitioner or therapist. The statements are strictly informative and may not be used as a substitute for a professional medi-cal treatment. Every user is beholden to ascertain by accurate examination and, if required, after consultation of a specialist, if the application is the right one for their specific cases. Each dose or application takes place at the risk of the user. We dissociate ourselves from any healing statements or promises.

Copyright pictures and text are subject to the copyright of Santé Sweden AB. and may not be used without explicit confirmation. Copyright © 2026 Santé Sweden AB. All rights reserved.

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Upphovsrättsskyddade bilder och text omfattas av Santé Sweden AB:s upphovsrätt och får inte användas utan uttrycklig bekräftelse. Copyright © 2026 Santé Sweden AB. Alla rättigheter förbehålls.